









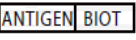

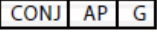



BioCLIA[®] Autoimmune Reagent Kit

GBM

Chemiluminescent Microparticle Immunoassay

Magnetic bead chemiluminescence immunoassay (CLIA) for semi-quantitative determination of anti-GBM antibody in human serum/plasma

Key to Symbols Used

	Catalog Number		Expiration Date
	For <i>In Vitro</i> Diagnostic Use		Lot Number
	Store between +2°C and +8°C		Consult Instruction for Use
	Manufacturer		Authorized Representative in European Union
	Biotinylated Antigen		Contains Sufficient for $< n >$ Tests
	Conjugate		Chemical Risk Warning
	Microparticle		Biological Risk Warning

BioCLIA® Autoimmune Reagent

Kit , GBM

Intended Use

BioCLIA GBM (Chemiluminescent Microparticle Immunoassay) is intended for the *in vitro* quantitative measurement of IgG antibodies directed to GBM in human serum as an aid in the diagnosis of vasculitis in conjunction with other laboratory and clinical findings. This kit is used on the instrument of BioCLIA® 1200 and BioCLIA® 6500.

Catalog Numbers

MY00113 (50 Tests/kit)

MY00164 (100 Tests/kit)

Summary and Explanation

Serological detection of anti-neutrophil cytoplasmic antibodies (ANCA) contributes to the autoimmune diseases diagnosis include Wegener's granulomatosis, acute progressive glomerulonephritis, polyarteritis, ulcerative colitis, and primary sclerosing cholangitis.^{1, 2} PR3, MPO and GBM are general indicators for the detection of ANCA, which can greatly improve the early diagnostic rate of renal vasculitis.

The main component of the glomerular basement membrane (GBM) is the extracellular matrix protein including type IV collagen, laminin, fibronectin and proteoglycans. The epitope of anti-GBM antibodies are located on the type IV collagen. Type IV collagen molecule is composed of three chains of 170 kDa. These chains form several triple-helix domains, and the domains are separated by the amino acid sequence which cannot form the helix. A tight spiral zone (7S domain) is located at the amino terminal and a spherical handle shaped structure (NC1 domain) at the carboxy terminal. The target antigen of anti-GBM antibodies is in NC1 domain of α -3 (IV) chain.

Anti-GBM antibody is a serological indicator for all anti-GBM glomerulonephritis including Goodpasture's syndrome.^{3, 4} In cases with no lung disorders, the positive rate of anti-GBM antibody is 60% while it is 80% - 90% in cases with lung disorders. Although the incidence of Goodpasture's syndrome is relatively low (only 0.5% of all kidney disorder patients), but the disease develop rapidly. If not treated well, the mortality rate will as high as 75-90%.^{5, 6} Early diagnosis and proper treatment can significantly reduce the mortality.

Principles of the Procedure

BioCLIA GBM assay is a two-step immunoassay using microparticle, enzyme-labeled chemiluminescent technology.

In the first step, the streptavidin coated magnetic

microparticle, the biotinylated GBM and human serum/plasma sample are mixed in an assay tube, which allows patient specific anti-GBM to bind. Secondly, after incubation, unbound biotinylated GBM and sample matrix are removed by washing, and the Microparticle-GBM-anti-GBM antibodies immune complexes are kept with the help of a magnetic separator. Third, anti-human IgG conjugated alkaline phosphatase is added. Fourth, after incubation, excess enzyme conjugates are removed by washing and finally the bound enzyme is detected by addition of chemiluminescent substrate. The relative light unit (RLU) intensity is proportional to the amount of anti-GBM specific IgG. According to the anti-GBM specific IgG RLU-concentration standard curve, the RLU tested can be interpreted to anti-GBM specific IgG concentration in the sample expressed as RU/mL.

For semi-quantitation of anti-GBM antibodies, the BioCLIA GBM assay utilizes a predefined lot specific Master Curve that is uploaded into the instrument through the reagent Master Calibration Curve barcode. Based on the Master Curve, and results obtained by running two Calibrators, an instrument specific Working Curve is created, which is used to calculate anti-GBM antibodies concentration RU/mL from the relative luminescent units (RLU) obtained for each sample.

Specimen Collection

The appropriate specimen types for BioCLIA GBM Reagents are human serum and plasma (Sodium citrate anticoagulant; Heparin anticoagulant; EDTA anticoagulant). Cloudy samples should be purified by low-speed centrifugation. To prevent erroneous results due to the presence of fibrin, ensure that complete clot formation has taken place prior to centrifugation of samples. Some samples, particularly those from patients receiving anticoagulant therapy, may require increasing clotting time.

Freshly collected specimens stored in refrigerator (2-8 °C) are valid for testing within 8 days. The stored specimen should reach to room temperature (18-25°C) before testing, and should not be stored in this temperature condition more than 2 days. All on board specimens should be tested within 10 hours. Three freeze (at -20°C) -thaw for specimens do not affect the testing results.

Warnings and Precautions

1. This assay is only for use in the BioCLIA® 1200 and BioCLIA® 6500.
2. This product requires the handling of calibrators, controls and human specimens which contain human sourced materials. It is recommended that all human sourced materials are considered to be potentially infectious and handled in accordance with the OSHA Standard on Bloodborne Pathogens.⁷ Biosafety Level 2 or other appropriate biosafety practices should be

used for materials that contain or possibly contain infectious agents.^{8, 9, 10} Avoid contacting with skin and eyes. Wear suitable protective gloves and clothing.

3. Liquid waste and solid waste are temporarily stored at separate containers. Waste management should also be handled in accordance with standards mentioned in chapter Warnings and Precautions point No. 2.

4. Spilled reagents should be cleaned up immediately. Observe all federal, state and local environmental regulations when disposing wastes.

5. Once opened, this reagent cartridge must be stored in the instrument's reagent carousel. Avoid spilling the reagents when the reagent cartridge is placed into the instrument.

6. Chemical contamination of the reagents can result from improper cleaning or rinsing of the instrument. Residues from common laboratory chemicals such as formalin, bleach, ethanol, or detergent can cause interference in the assay. Be sure to follow the recommended cleaning procedure of the instrument as outlined in the BioCLIA® 1200 and BioCLIA® 6500 generator.

7. Precautions: 

Proclin 300 added in the reagents (GBM Antigen, Conjugate) at concentration between 0.0015% - 0.6%.

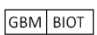
Storage Instructions

The kit is stable until the expiration date, if it is stored and handled. Routine store the kit in refrigerator (2-8 °C). Vial opened reagents or onboard reagents can be used up to 56 consecutive days (2-8 °C). The BioCLIA® 1200 and BioCLIA® 6500 software monitors the expiration of the reagent cartridge. The system will not accept expired reagents. Three freeze (at -20°C) -thaw cycles before testing has no effect on the kit reagents.

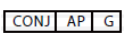
Materials Supplied

Components are matched in sets. Labels supplied within the kit will be needed for the assay testing.

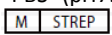
BioCLIA GBM

- **GBM Antigen** 1 bottle (2.5/5 mL) Biotinylated GBM antigen in 0.01 M PBS (pH7.4) buffer with stabilizer. 

Preservatives: 0.0015% < Proclin 300 < 0.6%

- **Conjugate** 1 bottle (6.75/13.5 mL) AP labeled anti-human IgG antibodies in 0.05 M MES (pH6.0) Buffer with stabilizer. 

Preservatives: 0.0015% < Proclin 300 < 0.6%.

- **Microparticle** 1 bottle (2.5/5 mL) Streptavidin-microparticles in 0.01 M PBS (pH7.4) buffer with stabilizer. 

Preservatives: 5-Bromo-5-Nitro-1, 3-Dioxane (BND) < 1%.

Kit Component Supplied Separately

Additional Materials Required But Not Provided:

- BioCLIA® 1200 (Cat No. MA00139)
- BioCLIA® 6500 (Cat No. MA00243)
- BioCLIA Autoimmune Calibrator Set, GBM (Cat No. MY00215, 2 x 1 mL; Cat No. MY00266, 4 x 1 mL)
- BioCLIA Autoimmune Control Set, GBM (Cat No. MY00317, 2 x 1 mL; Cat No. MY00368, 4 x 1 mL)
- BioCLIA Sample Diluent I (Cat No. MY00965)
- BioCLIA System Wash Buffer (Cat No. MY00404)
- BioCLIA System Substrate (Cat No. MY00405)
- BioCLIA Cuvettes (Cat No. MA00138, MA00244)
- BioCLIA Silicon Gasket (Small) (Cat No. MV00195)
- BioCLIA Silicon Gasket (Large) (Cat No. MV00196)
- BioCLIA Substrate Tube Maintenance cleanser (Cat No. MA00140)
- BioCLIA Sample Probe Maintenance cleanser (Cat No. MA00141)
- BioCLIA Micro Cup (Cat No. MA00142)
- Distilled Water

Assay Procedure

Note that, it is important to perform all routine maintenance procedures for optimal performance, such as routine cleaning, calibration and control procedures that are defined in the BioCLIA® 1200 and BioCLIA® 6500 User Manual.

See the BioCLIA® 1200 and BioCLIA® 6500 User Manual for preparation, setup, dilutions, adjustment, assay and quality control procedures.

Users should have the periodic calibration procedure for every 56 consecutive days from last calibration. Besides, a calibration procedure should be carried out when a new batch of BioCLIA GBM kit is used.

The control procedure should be done before running the specimens each day. Users also can adjust the control procedure period according to their own lab frequency.

Expected Values

Each Laboratory should establish its own reference ranges.

When the customer see a problem (High CV or unusual values, rerun controls and analyze specimens again.

Result Analysis

With the help of the built in master calibration curve and specified two-point calibrator set for the instrument, the BioCLIA will automatically calculate the auto-antibodies concentration of each specimen and interpret the results into RU/mL. The concentration of anti-GBM antibody sample is

reported as < 2 RU/mL when it is lower than the minimum detection limit, while reported as > 400 RU/mL when it is higher than the range of measurement.

Sample Dilution

The specimens are diluted with BioCLIA Sample Diluent I before testing (dilution ratio 1:20) by the BioCLIA® 1200 and BioCLIA® 6500 automatically.

Cut-Off Value Determination

120 clinical samples, including 30 positive sera, 30 negative sera, 30 positive plasma and 30 negative plasma were collected and valued. These samples were venous blood from human aged between 0 - 80, sealed and stored at 2 – 8 °C. Results of 120 clinical samples tested by the BioCLIA GBM kit were analyzed using the receiver-operating characteristic curve (ROC) and the cut-off value was determined at 20 RU/mL.

Test Result Interpretation

Specimen with concentration < 20 RU/mL, interpreted as negative;

Specimen with concentration ≥ 20 RU/mL, interpreted as positive.

Test results only reflect the sample collecting status and should be interpreted/analyzed for diagnosis in conjunction with other laboratory and clinical findings.

Performance Characteristics

APPEARANCE

Kit components are complete with no leakage. No precipitation or floc in liquid reagents. Packing labels are clear and easy to be identified.

ACCURACY / SPIKED RECOVERY

The accuracy/spiked recovery was determined by analyzing samples spiked with known amounts of anti-GBM antibodies into certain matrix. Anti-GBM antibody positive samples (low 100 RU/mL, mid 200 RU/mL, high 300 RU/mL) were spiked into two matrixes (50, 100 RU/mL) separately at the volume ratio of 1:9, making totally 6 spiked samples and each sample was tested in triplicate. The spiked recovery for the concentration of anti-GBM antibodies was calculated.*

Spike d Conc.	Matrix 50 RU/mL			Matrix 100 RU/mL		
	Obs	Exp	Obs/Exp	Obs	Exp.	Obs/Exp
Neat	49.5 4			103.3 5		
100 RU/mL	56.0 7	54. 6	102.7%	95.09	103. 0	92.3%
200 RU/mL	70.0 3	64. 6	108.4%	108.7 5	113. 0	96.2%
300 RU/mL	81.5	74.	109.3%	132.7	123.	107.9%

L	5	6		2	0	
---	---	---	--	---	---	--

*Representative data; results in individual laboratories may vary from these data.

TRACEABILITY

The reported values were determined with multiple runs on the BioCLIA® 1200 and BioCLIA® 6500 using specific reagents against an in-house standard. Results are reported in RU/mL, which is interpreted from relative light unit (RLU). Method comparison test showed good sensitivity and specificity.

PRECISION

A study based on guidance from CLSI document EP5-A2 was performed for determining the precision of BioCLIA GBM kit. Human serum in the in-house reference panel (RP1, RP2, RP3, RP4) was tested with 10 replicates per sample for intra-assay precision evaluation, while with 4 replicates per sample for inter-assay precision. Each sample tested in individual runs, and 2 runs per day for 10 days, a total of 80 points. Data from this study are summarized in the following table.*

Intra-assay precision: CV < 10%

Intra-Assay	RP1	RP2	RP3	RP4
Mean (RU/ml)	10.15	19.99	99.73	356.88
CV	1.8%	1.1%	1.2%	1.0%

Inter-assay precision: CV < 15%

Inter-Assay	RP1	RP2	RP3	RP4
Mean (RU/ml)	9.70	19.46	100.32	352.29
CV	4.9%	3.3%	3.5%	3.2%

*Representative data; results in individual laboratories may vary from these data.

LIMIT OF BLANK / DETECTION (LOB/LOD)

LOB/LOD was determined consistent with CLSI EP17-A guideline. LOB/LOD of the BioCLIA GBM assay was lower than 1.0 RU/mL, which is below the analytical measuring range of the assay.

ASSAY REPORTABLE RANGE

The BioCLIA GBM kit has a reportable linear range of 2 - 400 RU/mL. The linear range was determined by diluting a high positive anti-GBM antibody serum sample with a negative sample to several concentrations which covers the entire assay linear range according to the scheme in CLSI EP6-A. The expected value was plotted against the observed value, and linear regression analysis was performed to get slope, intercept and coefficient of correlation (r) values. The results are summarized in the tables below: *

Slope	Intercept	r
0.97	-0.28	0.99

Assay linear range is 2-400 RU/mL. Results below the lower limit will be reported as < 2 RU/mL, while those are above the upper limit will be reported as > 400 RU/mL.

*Representative data; results in individual laboratories may vary from these data.

INTERFERENCE

Bilirubin, hemoglobin, triglycerides, rheumatoid factor (RF) and human anti-mouse antibody (HAMA) will not affect the BioCLIA GBM assay performances when at the level indicated below.

Bilirubin ≤ 40 mg/dL;

Hemoglobin ≤ 150 mg/dL;

Triglycerides ≤ 1,000 mg/dL;

Rheumatoid factor (RF) ≤ 1,000 IU/mL;

Human anti-mouse antibody (HAMA) ≤ 2,000 ng/mL.

METHOD COMPARISON

Method comparison was implemented by comparing BioCLIA GBM assay to the predicated assay.

Clinical Sample		BioCLIA GBM		
		-	+	Total
Predicated Method	-	99	0	99
	+	0	1	1
	Total	99	1	100

Sensitivity	100.0%
Specificity	100.0%
Total agreement	100.0%

Limitations

- The effectiveness of this kit is only confirmed for human serum/plasma, the applicability of the other kinds of samples is not verified.
- Reliable and reproducible results will be obtained when the assay procedure is carried out in accordance with the instructions and with adherence to good laboratory practice.
- Clinical diagnosis should not be made on the findings of a single test result, but should be integrate with all clinical and laboratory findings.

References

1. Savige J, Davies D, Falk RJ, Jennette JC, Wiik A. Antineutrophil cytoplasmic antibodies and associated diseases: a review of the clinical and laboratory features. *Kidney International* 2000;57:846-62.
2. Bossuyt X. Serologic markers in inflammatory bowel disease. *Clinical Chemistry* 2006;52:171-81.
3. Hellmark T, Johansson C, Wieslander J. Characterization of anti-GBM antibodies involved in Goodpasture's syndrome. *Kidney International* 1994;46:823-9.
4. Salama AD, Levy JB, Lightstone L, Pusey CD. Goodpasture's disease. *Lancet* 2001;358:917-20.
5. Herody M, Bobrie G, Gouarin C, Grünfeld JP, Noel LH.

Anti-GBM disease: predictive value of clinical, histological and serological data. *Clinical Nephrology* 1993;40:249-55.

6. Kathuria P, Sanghera P, Stevenson FT, Sharma S, Lederer E, Lohr JW, Talavera F, et al. Goodpasture Syndrome Treatment & Management. Medscape Reference. 2013.

7. US Department of Labor, Occupational Safety and Health Administration, 29 CFR Part 1910.1030, Occupational Exposure to Bloodborne Pathogens. Jan 2001.

8. US Department of Health and Human Services. Biosafety in Microbiological and Biomedical Laboratories, Fourth Edition. Washington, DC: US Government Printing Office, May 1999.

9. World Health Organization. Laboratory Biosafety Manual. Geneva: World Health Organization.2004.

10. Clinical and Laboratory Standards Institute. Protection of Laboratory Workers from Occupationally Acquired Infections: Approved Guideline - Third Edition. CLSI Document M29-A3. Wayne, PA: Clinical and Laboratory Standards Institute, 2005.



HOB Biotech Group Co., Ltd



C6 Building, No. 218 Xinghu Road, Suzhou Industrial Park, Suzhou, Jiangsu, 215123, China

REGISTRANT/MANUFACTURE: HOB Biotech Group Co., Ltd

ADDRESS/LOCATION:

C6 Building, No. 218 Xinghu Road, Suzhou Industrial Park, Suzhou, Jiangsu, 215123 China

CONTACT INFORMATION: TEL (+86)512-69561996 Fax (+86)512-62956652

WEBSITE: www.hob-biotech.com

CUSTOMER SERVICE: HOB Biotech Group Co., Ltd

CUSTOMER SERVICE CONTACT: TEL (+86)4008601202



EUROPE REPRESENTATIVE: Emergo Europe

ADDRESS/LOCATION:

Prinsessegracht 20, 2514 AP The Hague, The Netherlands

Technical Assistance

For technical assistance, contact your National Distributor.

17th April 2019

Revision 9