

BioCLIA Autoimmune Calibrator Set, PCA

Cat.No.	Kit Size
MY00247	2 X 1 mL
MY00298	4 X 1 mL

INTENDED USE

The BioCLIA Autoimmune Calibrator Set, PCA is intended for the calibration of the BioCLIA PCA performed on the BioCLIA 6500 and BioCLIA 500.

For professional *in vitro* diagnostic use only.

SUMMARY AND EXPLANATION

Pernicious anemia (PA) is a disease in which there are not enough red blood cells partially due to autoimmune problems or lack of vitamin B₁₂. The most common initial symptom is feeling tired. Other symptoms may include shortness of breath, pale skin, chest pain, numbness in the hands and feet, poor balance, a smooth, red tongue, poor reflexes, and confusion. ¹ If treatment is not provided, some of these problems may become permanent. ²

PA may be considered as an end stage of immune gastritis, a disease characterized by stomach atrophy and the presence of antibodies to parietal cells and intrinsic factor. ³ This autoimmune disorder is localized to the body of the stomach, where parietal cells are located.

⁴ Antibodies to intrinsic factor and parietal cells cause the destruction of the oxyntic gastric mucosa, leading to the subsequent loss of intrinsic factor synthesis. Without intrinsic factor, the ileum can no longer absorb the B₁₂. ⁵

Parietal cell antibodies are found in other autoimmune disorders and also in up to 10% of healthy individuals, making the test nonspecific. However, around 85% of PA patients have parietal cell antibodies, which means they are a sensitive marker for the disease. ⁶ The combination of both tests of intrinsic factor antibodies and parietal cell antibodies may improve overall sensitivity and specificity of the diagnostic results. ⁷ About 90% of individuals with PA have antibodies for parietal cells; however, only 50% of all individuals in the general population with these antibodies have pernicious anemia. ⁸

MATERIALS SUPPLIED

- PCA Calibrator 1** Barcode labeled tubes with buffer containing human antibodies to PCA in stabilizers and preservatives. Ready to use, 1 mL.



Preservatives: 0.0015% < Proclin 300 < 0.6%.

- PCA Calibrator 2** Barcode labeled tubes with buffer containing human antibodies to PCA in stabilizers and preservatives. Ready to use, 1 mL.



Preservatives: 0.0015% < Proclin 300 < 0.6%.

The Calibrator Code contains calibrators' information is provided in each kit.

WARNINGS AND PRECAUTIONS

- For professional *in vitro* diagnostic use only.

- Do not use any calibrators beyond their expiration dates.
- Do not mix calibrators from different lots unless specified.
- Instructions must be carefully followed for using and storing of calibrators. Any modification in procedure may interfere with the results.
- Calibrators and contaminated vials must be handled strictly following safety guidelines or rules of biological hazards to ensure the users' and environmental safety.
- Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the competent authority of the Member State in which the user and/or the patient is established!

Precautions:



Human serum is added in the calibrators.

- The human derived material in this product was tested by FDA approved methods and found nonreactive for Hepatitis B Surface Antigen (HBsAg), Anti-HCV and HIV 1/2 antibodies. Handle as if potentially infectious. ⁷ Avoid contacting with skin and eyes. Do not empty into drains. Wear suitable protective gloves and clothing.



- Proclin 300 is added in the calibrators at concentration between 0.0015% - 0.6%.
- Calibrators contain chemical and biological components. Avoid ingesting or splashing onto skin and mucous membrane. If direct contact with calibrators happens, rinse immediately the contact surface with plenty of water and see a doctor if necessary.

STORAGE CONDITIONS

- Store the kit at 2-8 °C.
- The shelf life of the unopened kit is 12 months from day of production.
- Vial opened calibrators could be used for 28 successive days, exposure no more than 2 hours each time when kept uncapped and is good for up to 20 calibrations, after which the reagent must be discarded.
- Avoid repeated freezing and thawing.

ASSAY PROCEDURE

Detailed information about operating the BioCLIA instruments can be

taken from the Instrument User's Manual.

Note that, it is important to perform all routine maintenance procedures for optimal performance.

Assay Calibration

The BioCLIA Autoimmune Reagent Kit utilizes a predefined lot specific Master Curve which is uploaded into the instrument via the barcode provided in the main reagent kit. The Calibrator Code contains calibrator information is then scanned. Based on the results of running two calibrators, the instrument specific Working Calibration Curve is generated and is used to calculate the concentration from the RLU obtained for each patient.

For each new lot of reagent, please calibrate prior to the first time use, and every 28 days thereafter. The software will not allow the lot to be used if the above requirements are not meet.

Programming and Running samples

1. Put the kit into the corresponding position of the reagent chamber of the fully automatic chemiluminescence analyzer. The information of the kit can be uploaded into the instrument system through the scanning of reagent barcode, and can also be set through the supporting software of the instrument.
2. The information of calibrator / quality control is identified by scanning the calibrator / control barcodes, and the position of calibrator / quality control is assigned in the instrument system.
3. The sample to be tested is placed on the instrument sample rack chamber, and the corresponding test information is edited through the instrument supporting software.
4. Start the operation procedure, and all calibrator / quality control / sample processing steps will be automatically executed.

TRACEABILITY

The reported values were determined over multiple runs on the BioCLIA 6500 and BioCLIA 500 using specific lots of reagents against an in-house standard. DGP IgA results are reported in RU/mL which is interpreted from relative light unit (RLU). Method comparison test showed good sensitivity and specificity of tested assay.

LIMITATIONS

- The calibrators are designed for calibration of the same lot of BioCLIA Autoimmune Reagent Kit.
- The calibrators can be kept uncapped onboard the instrument up to 2 hours for each time of usage. And a total up to 20 calibrations are suggested, for any longer period of time, the reagent should be discarded, otherwise may result in improper calibration of the assay and which can give improper results.

SYMBOLS

	Calibrator 1
	Calibrator 2

	Catalog Number		Use-by date
	In Vitro diagnostic medical device		Lot Number
	Store between +2°C and +8°C		Consult Instruction for Use
	Manufacturer		Authorized Representative in the European Community
	CE Marking		Contains Sufficient for <n> Tests
	Biological Risk		GH507 Warning

REFERENCE

1. Toh BH, van Driel IR, Gleeson PA. Pernicious anemia. New England Journal of Medicine 1997;337:1441-8.
2. Hsing AW, Hansson LE, McLaughlin JK, Nyren O, Blot WJ, Ekblom A, Jr FJ. Pernicious anemia and subsequent cancer. A population-based cohort study. Cancer 1993;71:745-50.
3. De Block CE, De Leeuw IH, Van Gaal LF. Autoimmune gastritis in type 1 diabetes: a clinically oriented review. Journal of Clinical Endocrinology & Metabolism 2008;93:363-71.
4. Banka S, Ryan K, Thomson W, Newman WG. Pernicious anemia - genetic insights. Autoimmunity Reviews 2011;10:455-9.
5. Lahner E, Annibale B. Pernicious anemia: New insights from a gastroenterological point of view. World Journal of Gastroenterology 2009;15:5121-28.
6. Snow CF. Laboratory diagnosis of vitamin B12 and folate deficiency: a guide for the primary care physician. Archives of Internal Medicine 1999;159:1289-98.
7. Gräsbeck R. Imerslund-Gräsbeck syndrome (selective vitamin B12 malabsorption with proteinuria). Orphanet Journal of Rare Diseases 2006;1:4677-87.
8. Butler CC, Vidal-Alaball J, Cannings-John R, McCaddon A, Hood K, Papaioannou A, et al. Oral vitamin B12 versus intramuscular vitamin B12 for vitamin B12 deficiency: a systematic review of randomized controlled trials. Family practice 2006;23:279-85.
9. Richmond JY, McKinney RW. Biosafety in microbiological and biomedical laboratories: U.S.GPO.1999.



HOB Biotech Group Corp., Ltd.
C6 Building, No. 218 Xinghu Road, Suzhou Industrial Park,
Suzhou, Jiangsu, 215123, China

CONTACT INFORMATION:


TEL (+86)512-69561996

Fax (+86)512-62956652

WEBSITE: www.hob-biotech.com

CUSTOMER SERVICE: HOB Biotech Group Co., Ltd

CUSTOMER SERVICE: TEL (+86)4008601202

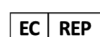

EUROPE REPRESENTATIVE: Emergo Europe

ADDRESS/LOCATION :

Prinsessegracht 20, 2514 AP The Hague, The Netherlands



The eIFU is available on Website:

<http://en.hob-biotech.com/usercenter/login.aspx>

TECHNICAL ASSISTANCE

For technical assistance, contact your National Distributor.

Date of issue: 17th March 2019

Date of revision: 20th December 2021

Change Control Number: CN21129E

Version: A/1 (EN)